

SPECIFICATION

SERVOspin Next Table Top Centrifuge



Easy-to-learn operating concept with programme memory function

The SERVOspin Next has an easy-to-learn operating concept. The speed can be set manually in steps of 100 rpm, and the running time can be selected in minutes up to a maximum of 99 minutes. To simplify operation, 5 freely selectable programmes are available

Centrifugation of sample tubes

The SERVOspin Next benchtop centrifuge is ideally suited for centrifuging test tubes containing a separation medium such as gel or granulate. Due to the swing-out rotor, centrifugation takes place in the direction of the tube axis, so that the centrifuged components form horizontally uniformly thick layers.




Centrifugal force with clear operating concept

The clear operating concept is easy to learn and offers the user the greatest possible comfort and safety. For the daily routine, 5 freely programmable programme slots are available. The powerful motor with fan cooling and the integrated unbalance control ensure very smooth and quiet operation. For the user's safety, the lid is electromagnetically locked. During centrifugation, the speed and the remaining time are shown on the display. The display can be switched to show the RCF (centrifugal force) by simply pressing a button.

The advantages at a glance

- Versatile use in laboratory, doctor's office, hospital, pharmacy, etc.
- Oscillating rotor 6 x 15 ml with adapter set for adaptation to different tubes.
- Ergonomic control panel and display.
- Digital preselection of running time and speed.
- Display switchable from speed (rpm) to centrifugal force (rcf).
- Electromagnetic lid lock.
- Emergency opening in case of power failure .
- Sturdy housing and long-life motor.
- Easy rotor change.
- Automatic unbalance control.
- Easy cleaning and disinfection.
- 24 months warranty.

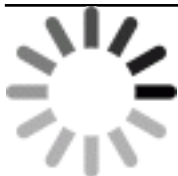

Technical data

- Speed: 100 - 4000 rpm
 - Centrifugal force (Max.RCF): 2510 G
 - Running time: 1 to 99 min
 - Dimensions: 420 X 340 X 280 mm (D x W x H)
 - Weight: 12 kg
 - Mains voltage: 220 - 230 Volt 50/60 Hz
-

SPECIFICATION

SERVOspin Next Table Top Centrifuge

Type	PZN	REF	Unit
SERVOspin Next incl. swing rotor		H7 820	1 piece

Product	REF	Quantity	Unit price	Total

_____ plus tax

Personal details:

Customer number: _____

Company: _____

First name: Last name: _____

E-mail: _____

Street, street#: _____

Postal code, city: _____

I hereby place a binding order:

Company stamp / Address

_____ Date / Sign

You are ordering under our terms and conditions, our valid terms and conditions can be found at <https://www.servoprax.de/en/terms-and-conditions/>